


Feroza Dengue NS1 Ag Test

One step qualitative immunochromatographic test for detection of dengue virus NS1 antigen in human serum, plasma or whole blood

ISO 9001:2008, ISO 13485:2003, Thai GMP Certified

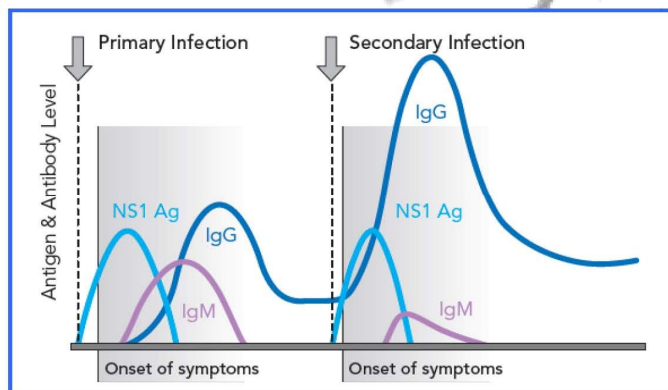

- Only one step method
- Store at room temperature without refrigeration

Test Procedure

Take a device from the pouch and place it on a flat.

Add 3 drops or 100 µl of serum, plasma or whole blood into the sample well.

Interpret the test result within 15-20 minutes

Interpretation

Positive	
Negative	
Invalid	

Reference table of result and its interpretation

Dengue NS1 Ag	Dengue Results		Interpretation
	Dengue IgM	Dengue IgG	
-	-	-	Negative for dengue infection
+	-	-	Acute primary dengue, early phase
+	+	-	Acute early dengue infection
-	+	-	Primary dengue infection
+	-	+	Secondary acute infection
+	+	+	Repeated dengue infection
-	+	+	Acute dengue infection, repeated infection
-	-	+	Passed infection or late of secondary infection

Note : For the most accurate dengue diagnosis in patients presenting at various stages of illness, a combination of the Feroza Dengue IgM/IgG Duo Test and Feroza Dengue NS1 Ag Test should be used.

Ordering information :-

Product code : DEN-02-025 Feroza Dengue NS1 Ag Test(25 Tests/Kit)

Learn more about our products on our web site : www.brianet.com



Manufactured by : Pacific Biotech Co., Ltd.
42 Moo 4, Petchaboon-Chianglue Rd., Muang, Petchaboon, Thailand 67000
Tel : 6656-748-650-1 Fax : 6656-748-649
<http://www.brianet.com> ; email: tosapol@brianet.com

Distributed by : Pacific Biotech Co., Ltd.
6 Sonthiwatana 3, Ladprao 110, Wangthonglang, Bangkok, Thailand 10310
Tel : 662-5302760 ; 662-9351281 Fax: 662-9318340
<http://www.brianet.com> ; email: tosapol@brianet.com